



## AMERICAN INTERNATIONAL UNIVERSITY - BANGLADESH

Department of Natural Science (Chemistry)

Faculty of Science &amp; Technology

Programs: B.Sc. Eng'g (EEE/CSE/IPE)

CHEM 1101: CHEMISTRY

## Chemistry Lab Report

Semester: Summer

Session: 2021-2022

NO EXPERIMENT, NO REPORT

Experiment No: 05

Name of the Experiment: Standardization of Sodium thiosulphate ( $\text{Na}_2\text{S}_2\text{O}_3$ ) solution with standard potassium dichromate ( $\text{K}_2\text{Cr}_2\text{O}_7$ ) solution

Date of Performance: 05-07-22, Date of Submission: 19-07-22

Course-Teacher: Dr. Saiful Islam

## Instructions:

1. A lab report consists of three parts: a cover page, body of the report and a data and results sheet (lab-sheet).
2. This is the cover page of a report and students will collect and preserve the lab-sheet of a particular experiment to be performed.
3. Body of the report includes-(1) Objective of the Experiment, (2) Theory, (3) Name of the Chemicals, (4) Name of the Apparatus, (5) Percentage of Error (if necessary) and (6) Discussion (I. Precautions taken, II. Possible errors).
4. Use A1-size off-set paper, write on one side of the paper by hand keeping suitable margin.
5. Staple the lab-sheet at the end of the report and cover page on the top.
6. Submit the report in time to avoid deduction of marks.
7. Students working in a group will write and submit the report individually.
8. Copying of the report from others is strictly prohibited.

Name of the Student: Khondoker Md. Sabir Hasan

ID No: 21-45306-2, Section: M, Group: 9

FOR FACULTY USE ONLY

Faculty comments: \_\_\_\_\_, Signature: \_\_\_\_\_, Date: \_\_\_\_\_

Objective: To know the strength of  $\text{Na}_2\text{S}_2\text{O}_3$  solution (being a solution made from secondary standard substance) against standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution by oxidation-reduction titration.



Theory:

(i) Method:



Redox titration (also called oxidation-reduction titration) is a type of titration based on a redox reaction between the analyte (substance subjected to analysis) and titrant (standard solution taken in burette). Redox titration may involve the use of a redox indicator and/or a potentiometer.

(ii) Reaction:

$\text{Na}_2\text{S}_2\text{O}_3 \cdot 5\text{H}_2\text{O}$  is a secondary standard substance, a reducing agent, white solid and colorless in water solution.  $\text{K}_2\text{Cr}_2\text{O}_7$  is a primary standard substance, an oxidizing agent, orange solid and orange color in water solution. The balanced reactions of  $\text{K}_2\text{Cr}_2\text{O}_7$  in presence of  $\text{HCl}$  and  $\text{KI}$  (white solid, a reducing agent) and of  $\text{Na}_2\text{S}_2\text{O}_3$  and  $\text{I}_2$  are as follows:

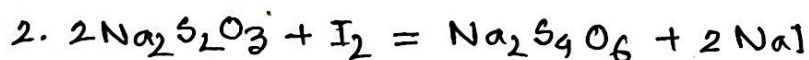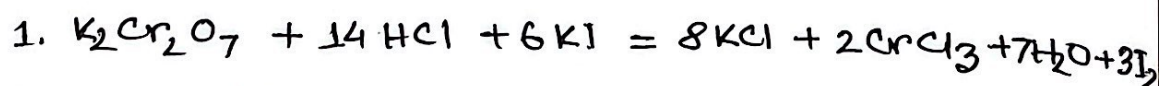

Redox Half Reactions,

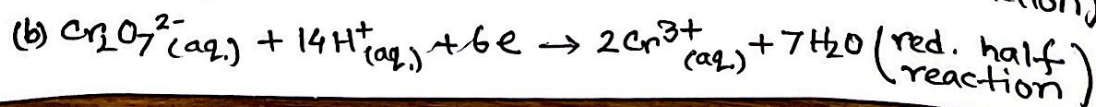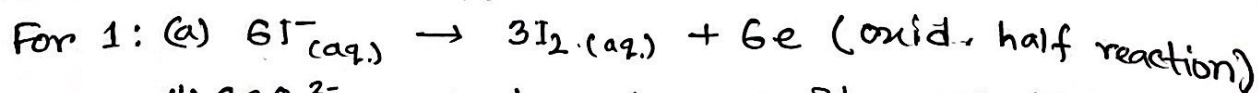



For 2: (a)  $2S_2O_3^{2-} (aq.) \rightarrow S_4O_6^{2-} + 2e^-$  (oxid. half reaction)

(b)  $I_2 (aq.) + 2e^- \rightarrow 2I^- (aq.)$  (red. half reaction)

### (iii) Indicator:

Starch solution is used in this titration involving iodine because it forms an intense blue complex with even a trace of iodine. But starch is not a redox indicator; it responds specifically to the presence of  $I_2$ , not to a change in redox potential. The active fraction of starch is amylose, a polymer of the sugar  $\alpha$ -D-glucose. In the presence of starch, iodine forms  $I_5^-$  chains inside the amylose helix and the color turns dark blue.

### Apparatus:

Burette (50mL), pipette (10mL), conical flask (250mL), volumetric flask (100mL), watch glass, pipette filler, dropper, stand and clamp etc.

### Required chemicals:

1. Potassium iodide.  $[KI]$
2. Sodium bicarbonate.  $[NaHCO_3]$
3. Conc. Hydrochloric acid.  $[HCl]$
4. Potassium dichromate.  $[K_2Cr_2O_7]$
5. Sodium thiosulphate.  $[Na_2S_2O_3]$
6. Starch.  $[(C_6H_{10}O_5)_n]$

Name: Khondoker Md. Sabit Hasan ID No: 21-45306-2 Section (Group): M (9)

AMERICAN  
INTERNATIONAL  
UNIVERSITY –  
BANGLADESH (AIUB)



# Experiment 5

CHEM 1101: CHEMISTRY (EEE/CoE/CSE/IPE)

## EXPERIMENT NO. 5: STANDARDIZATION OF SODIUM THIOSULPHATE ( $\text{Na}_2\text{S}_2\text{O}_3$ ) SOLUTION WITH STANDARD POTASSIUM DICHROMATE ( $\text{K}_2\text{Cr}_2\text{O}_7$ ) SOLUTION.

**OBJECTIVE:** To know the strength of  $\text{Na}_2\text{S}_2\text{O}_3$  solution (being a solution made from secondary standard substance) against standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution by oxidation-reduction titration.

### THEORY:

- (i) Method: Redox titration
- (ii) Reaction:
  1.  $\text{K}_2\text{Cr}_2\text{O}_7 + 14\text{HCl} + 6\text{KI} = 8\text{KCl} + 2\text{CrCl}_3 + 7\text{H}_2\text{O} + 3\text{I}_2$
  2.  $2\text{Na}_2\text{S}_2\text{O}_3 + \text{I}_2 = \text{Na}_2\text{S}_4\text{O}_6 + 2\text{NaI}$
- (iii) Indicator: Starch solution

### APPARATUS:

Burette (50mL), pipette (10mL), conical flask (250mL), volumetric flask (100mL), watch glass, pipette filler, dropper, Stand and clamp etc.

### REQUIRED CHEMICALS:

1. 12% KI solution,
2.  $\text{NaHCO}_3$ ,
3. Conc. HCl acid,
4. Standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution
5.  $\text{Na}_2\text{S}_2\text{O}_3$  solution,
6. Starch solution

### PREPARATION OF APPROX. 0.1N POTASSIUM DICHROMATE SOLUTION.

Transfer approx. 0.49 gram of pure  $\text{K}_2\text{Cr}_2\text{O}_7$  into a 100 mL measuring flask and then dissolve it with distilled water up to the mark.

$$\begin{aligned} \text{Strength of } \text{K}_2\text{Cr}_2\text{O}_7 \text{ solution} &= \frac{\text{Weight taken (in gm)} \times 0.1}{0.49} \text{ (N)} \\ &= \frac{0.48 \times 0.1}{0.49} \text{ (N)} \\ &= 0.098 \text{ (N)} \end{aligned}$$

May, 2022

Chemistry Lab Sheet



### PROCEDURE:

(Expt. 5 contd.)

Take 4 mL of 12% KI solution in a conical flask and dilute to about 50 mL. Add about one gm of  $\text{NaHCO}_3$  and shake the flask until the salt dissolves. Add 4 mL conc. HCl acid and then add 10 ml standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution by means of a pipette in the same flask. Shake the flask and cover it with a watch glass, allow the solution to stand for about five minutes in the dark (inside the desk or dark chamber). Rinse the watch glass and dilute the solution about 100mL. Titrate the liberated iodine with sodium thiosulphate solution from a burette until the brown color fades (light yellow). Add about 1 mL starch solution and continue titration by adding sodium thiosulphate solution from the burette until one drop of the sodium thiosulphate solution changes the color of the solution from deep blue to light green. This is the end point. Repeat the whole experiment 2-3 times. Calculate the strength of sodium thiosulphate solution.

### EXPERIMENTAL DATA:

Table: Standardization of supplied  $\text{Na}_2\text{S}_2\text{O}_3$  solution against standard  $\text{K}_2\text{Cr}_2\text{O}_7$  solution by oxidation-reduction titration.

| No. of reading | Vol. of $\text{K}_2\text{Cr}_2\text{O}_7$ (in mL) | Vol. of $\text{Na}_2\text{S}_2\text{O}_3$ (burette reading) (in mL) |       |            | Mean (in mL) |
|----------------|---|---|-------|------------|--------------|
|                |   | Initial   | Final | Difference |              |
| 1              | 10  | 0   | 7.8   | 7.8        | 7.6          |
| 2              | 10  | 7.8   | 15.1  | 7.3        |              |
| 3              | 10  | 15.1  | 22.8  | 7.7        |              |

### CALCULATIONS:

Strength of supplied  $\text{Na}_2\text{S}_2\text{O}_3$  solution:

$$V_{\text{Na}_2\text{S}_2\text{O}_3} \times N_{\text{Na}_2\text{S}_2\text{O}_3} = V_{\text{K}_2\text{Cr}_2\text{O}_7} \times N_{\text{K}_2\text{Cr}_2\text{O}_7}$$

$$\Rightarrow 7.6 \text{ ml} \times N_{\text{Na}_2\text{S}_2\text{O}_3} = 10 \text{ ml} \times 0.098 \text{ N}$$

$$\therefore N_{\text{Na}_2\text{S}_2\text{O}_3} = 0.129 \text{ N}$$

### RESULTS:

$\therefore$  Strength of the supplied  $\text{Na}_2\text{S}_2\text{O}_3$  solution is 0.129 N.

Students should know

- What are redox reaction, oxidizing agent and reducing agent?
- What is the difference between acid-base and redox indicator?
- Why it is necessary to keep your experimental solution in the dark?
- Is it iodometric or iodimetric that you are performing?
- Tell molecular weight and gram equivalent weight of  $\text{K}_2\text{Cr}_2\text{O}_7$  and  $\text{Na}_2\text{S}_2\text{O}_3$ .
- What is the function of starch?
- Can you calculate the normality and molarity of  $\text{K}_2\text{Cr}_2\text{O}_7$  and  $\text{Na}_2\text{S}_2\text{O}_3$ ?

Text: M. Mahbubul Huque and A. Jabber Mian, "Practical Chemistry", 2<sup>nd</sup> ed. (1972)

May, 2022

Chemistry Lab Sheet

## Discussion:

### (a) Precautions taken:

1. Avoid skin and eye contact with the chemicals.
2. Identify the safety equipment.
3. Wear clothing that covers torso and legs.

### (b) Possible errors:

~~I might added more than 100ml distilled water in the HCl acid solution~~

1. Weight was not taken properly.
2. Error might be occurred while taking the burette reading.
3. Using the equipment incorrectly.